



Lab Accessories

PHA-EH and PHS-11 Series
Starts at
\$69

PHA-EH Electrode Holder

- ✓ Hands Free Measurement
- ✓ Weighted Base for Stability
- ✓ Holds Multiple Sensors up to 14 mm (0.55") Diameter

PHS-11 Stirrer

- ✓ Variable Speed
- ✓ 120 Vac Adaptor Included




PHA-EH, \$69, shown smaller than actual size.

The PHS-11 Series of magnetic stirrers operate by generating a rotating electromagnetic field in the stirring heads. The liquid to be stirred has a magnetized stirring bar placed in it. The bar and liquid is rotated by the magnetic field produced in the heads. The PHS-11C also includes an electrode holder for hands-free measurements. The PHA-EH is an electrode holder only for applications where stirring is not desired. It has a weighted base for stability and locking knobs to hold the desired position.


Specifications

Working Dia: 120 mm (4 3/4")
RPM Range: 0 to 2300
Power: 110 or 220 Vac
Weight: 460 g (1lb)
Stirrer: PTFE coated magnet
Rod Length: 250 mm (10")
Base: 115 x 150 mm (4 1/2 x 6")

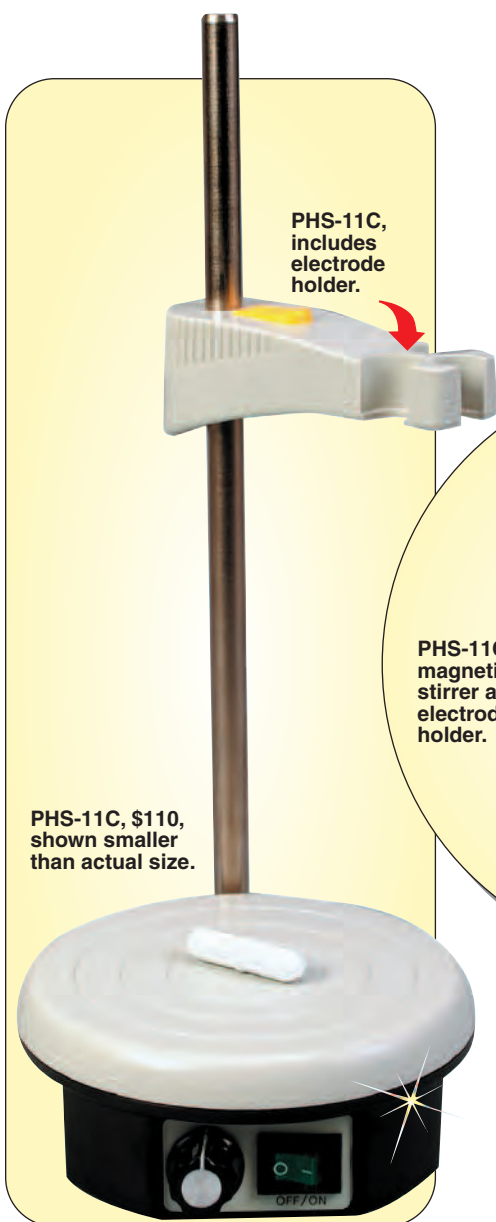

PHS-11C, includes electrode holder.

PHS-11C, \$110, shown smaller than actual size.

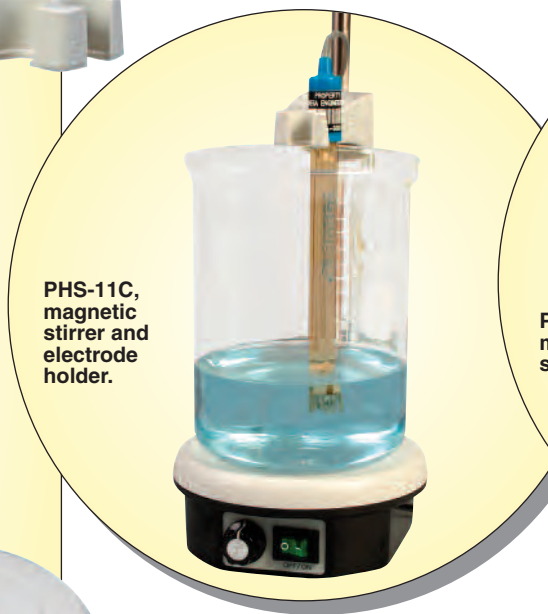

PHS-11C, magnetic stirrer and electrode holder.



PHS-11 magnetic stirrer.

AVAILABLE FOR FAST DELIVERY!

To Order (Specify Model Number)

Model No.	Price	Description
PHA-EH	\$69	pH electrode holder
PHS-11	105	Magnetic stirrer
PHS-11C	110	Magnetic stirrer with electrode holder

Comes complete with stirrers, power supply and operator's manual.

Ordering Example: PHS-11C, magnetic stirrer with electrode holder, \$110.

Order Online

omega.com

Over 100,000 Products Available!

To download information and to order pH Field and Lab Instrument products online, visit omega.com



UNITED STATES

www.omega.com
1-800-TC-OMEGA
Stamford, CT.

CANADA

www.omega.ca
Laval(Quebec)
1-800-TC-OMEGA

GERMANY

www.omega.de
Deckenpfronn, Germany
0800-8266342

UNITED KINGDOM

www.omega.co.uk
Manchester, England
0800-488-488

FRANCE

www.omega.fr
Guyancourt, France
088-466-342

CZECH REPUBLIC

www.omegaeng.cz
Karviná, Czech Republic
596-311-899

BENELUX

www.omega.nl
Amstelveen, NL
0800-099-33-44



More than 100,000 Products Available!

• Temperature

Calibrators, Connectors, General Test and Measurement Instruments, Glass Bulb Thermometers, Handheld Instruments for Temperature Measurement, Ice Point References, Indicating Labels, Crayons, Cements and Lacquers, Infrared Temperature Measurement Instruments, Recorders Relative Humidity Measurement Instruments, RTD Probes, Elements and Assemblies, Temperature & Process Meters, Timers and Counters, Temperature and Process Controllers and Power Switching Devices, Thermistor Elements, Probes and Assemblies, Thermocouples Thermowells and Head and Well Assemblies, Transmitters, Wire

• Flow and Level

Air Velocity Indicators, Doppler Flowmeters, Level Measurement, Magnetic Flowmeters, Mass Flowmeters, Pitot Tubes, Pumps, Rotameters, Turbine and Paddle Wheel Flowmeters, Ultrasonic Flowmeters, Valves, Variable Area Flowmeters, Vortex Shedding Flowmeters

• pH and Conductivity

Conductivity Instrumentation, Dissolved Oxygen Instrumentation, Environmental Instrumentation, pH Electrodes and Instruments, Water and Soil Analysis Instrumentation

• Data Acquisition

Auto-Dialers and Alarm Monitoring Systems, Communication Products and Converters, Data Acquisition and Analysis Software, Data Loggers Plug-in Cards, Signal Conditioners, USB, RS232, RS485 and Parallel Port Data Acquisition Systems, Wireless Transmitters and Receivers

• Pressure, Strain and Force

Displacement Transducers, Dynamic Measurement Force Sensors, Instrumentation for Pressure and Strain Measurements, Load Cells, Pressure Gauges, Pressure Reference Section, Pressure Switches, Pressure Transducers, Proximity Transducers, Regulators, Strain Gages, Torque Transducers, Valves

• Heaters

Band Heaters, Cartridge Heaters, Circulation Heaters, Comfort Heaters, Controllers, Meters and Switching Devices, Flexible Heaters, General Test and Measurement Instruments, Heater Hook-up Wire, Heating Cable Systems, Immersion Heaters, Process Air and Duct, Heaters, Radiant Heaters, Strip Heaters, Tubular Heaters